



## Introducing the SureTemp<sup>®</sup> Plus 692 Electronic Thermometer

Speed and accuracy you require.  
Safety and security you demand.

**WelchAllyn<sup>®</sup>**

Advancing Frontline Care<sup>™</sup>



## The SureTemp® Plus 692 – with the speed and accuracy you require, and the safety and security you demand

- **Easy to use, intuitive operation**
- **Removable probe well reduces the risk of cross-contamination**
- **Improved accuracy for even more reliable readings**
- **Innovative security features that cut down on the chance of theft**

When we set out to improve the Welch Allyn SureTemp, we knew we had our work cut out for us. After all, when you start with the thermometer that helped set the standard for speed and accuracy, the expectations run pretty high.

So how did we do? See for yourself by taking a look at the innovative SureTemp Plus 692 – easier, safer and more accurate than the original, and still built to withstand constant use in demanding environments.

**Easy:** An easy-to-read LCD display with flashing icons that clearly communicate in any language, and an intuitive user interface and ergonomic shape make using the SureTemp Plus 692 as simple as possible.

**Safe:** Patient safety is improved thanks to an interchangeable probe well that reduces the risk of cross-contamination.

**Fast and Accurate:** In addition to improved accuracy, the SureTemp 692 is one of the fastest thermometers on the market, providing approximately 4- to 6-second oral, 10-second paediatric axillary, 15-second adult axillary, and 10-second rectal measurements.

**Secure:** To top it all off, our new security features make sure the unit stays where it is supposed to be. A location ID on the display and other innovative features help prevent theft.

ID Field for electronic labelling of thermometer



Battery life indicator

User selectable icon and/or words for identification of thermometer mode

Interchangeable probe with push-button probe cover ejection

Removable probe well for easy cleaning

Convenient storage for 25 probe covers

C/F button for toggling between Fahrenheit and Celsius

Large, easy-to-read, backlit LCD

Mode button for selecting oral, paediatric axillary, or adult axillary

Button for timing pulse and respiration

Ergonomic design with rubberised grips


Prominent recall button for display of last temperature taken

Storage space for 50 additional probe covers in the wall mount (Wall mount standard with the 692 model)



Optional rolling stand for security and portability

Interchangeable rectal probe and well reduces risk of cross-contamination

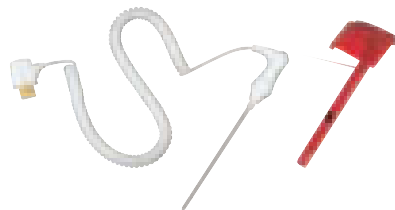

# The SureTemp® Plus 692

**WelchAllyn®**

Advancing Frontline Care™

The SureTemp Plus 692 –  
 setting the standard of  
 care everywhere

## SureTemp Plus 692 specifications

Calibration Accuracy	±0.1°C (Monitor Mode)
Temperature Predict Time	
Oral	approx. 6 secs
Adult Oesophageal (16 years and older)	approx. 2.5 secs
Paediatric Oesophageal (16 years and younger)	approx. 30 secs
Rectal	approx. 30 secs
Unit Dimensions	25 mm x 45 mm x 62 mm
Weight	15 grams
Ambient Temperature Range	0°C to 40°C
Patient Temperature Range	26°C to 42°C
Power Source	(3) 1.5V dc batteries
Battery Life	approx. 5000 temperatures
Warranty	1 year limited instrument, 1 year limited probe, and 90 day limited probe well



## Ordering information

### Thermometers

- SureTemp Plus Electronic Thermometer (Model 692), Wall Mount, Security System with ID Location Field, and Oral Probe with Oral Probe Well, 1 year Limited Warranty
- SureTemp Plus Electronic Thermometer (Model 692), Wall Mount, Security System with ID Location Field, and Rectal Probe with Rectal Probe Well, 1 year Limited Warranty
- SureTemp Plus Electronic Thermometer (Model 692), One Per Room, Wall Mount, 25 mm Oral Probe with Oral Probe Well, 1 year Limited Warranty *to be located to the wall*
- SureTemp Plus Electronic Thermometer (Model 692), One Per Room, Wall Mount, 25 mm Rectal Probe with Rectal Probe Well, 1 year Limited Warranty *to be located to the wall*

### Accessories

- Wall holder
- Wall holder for Spare Probe Well
- Calibration Key Assy, 690/692
- Probe Well Kit 25 mm Rectal
- Probe Well Kit 25 mm Rectal
- Probe Well Kit 25 mm Oral
- Probe Well Kit 25 mm Oral
- Probe Well, Blue
- Probe Well, Red
- Thermometer Stand with Wall holder for Spare Probe Well
- Bracket, M600 Stand Upgrade, with Wall holder for Spare Probe Well
- Extended Warranty M692

### Replacement probe covers

- Case of 1,000 probe covers

Welch Allyn UK Ltd  
 Winton Road, Stonehouse, Gloucestershire  
 Gloucestershire PL22 7YU, England  
 Tel: 0206 669 2696

Welch Allyn European Customer Service  
 Winton Road, Stonehouse, Gloucestershire  
 Gloucestershire PL22 7YU, Ireland  
 Tel: 0696 92696

Welch Allyn PT LTD, South Africa  
 9 Surrey Road, Durban 29, South Africa  
 P.O. Box 9, Montclair, South Africa 22  
 Tel: 031 269 696

Call Welch Allyn to order the SureTemp® Plus 692  
 State Street Road, P.O. Box 22,  
 San Mateo, CA 94062, USA  
 Tel: 666 666 666  
[www.welchallyn.com](http://www.welchallyn.com)

**WelchAllyn®**

Advancing Frontline Care™