



Hillrom™



The Hillrom TS7000dV Surgical Table is designed to work with Integrated Table Motion for the *da Vinci Xi*® Surgical System from Intuitive Surgical.

# TS7000DV SURGICAL TABLE

## REDEFINING PATIENT POSITIONING IN *DA VINCI*® SURGERY

Integrated Table Motion (ITM) enables the *da Vinci Xi*® System to connect to Hillrom's TS7000dV Surgical Table so the patient can be dynamically positioned while the surgeon operates. This allows the surgical team to efficiently manage the patient and utilize gravity for optimal exposure and access to the target anatomy.



## TS7000dV OR Table Package

The TS7000dV Surgical Table is the central piece of a unique package which provides components and accessories for multi-specialty *da Vinci*® surgery:

- Double Joint Head Section
- Two Parts Leg Section, Spreadable
- dV Remote Control
- dV Connection Cable

With any TS7000dV Operating Table two years of comprehensive service coverage and maintenance is included.



## Dramatic Range of Motion


Patient Weight Capacity	450 kg (1,000 lbs)	Exceptional weight capacity in any articulation, including slide
Height Range <sup>1</sup>	59 - 114 cm (±1 cm) (23 - 45"(± 0.4"))	Ultra low and high platform for the surgical team, including the surgical robot
Trendelenburg / reverse Trendelenburg	± 45°	Superior positioning to expose target anatomy, including during robotic-assisted surgery

## Additional Advanced Positioning Accessories

Hillrom offers an additional accessory package of specialty components that optimize patient positioning capabilities for complex robotic-assisted surgeries.



Two-Part Leg Section for midline access.

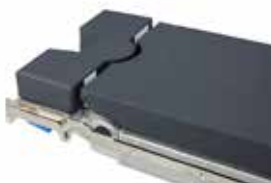

Pelvic Extension for enhanced imaging capability.



TrenGuard® for secure, steep Trendelenburg positioning.

For more information please contact your local Hillrom representative.

**hillrom.com**

Hillrom  
Amerika Building (7th Floor), Hoogoorddreef 15  
1101 BA Amsterdam, The Netherlands

da Vinci and da Vinci Xi are trademarks or registered trademarks of Intuitive Surgical, Inc. Trenguard is a registered trademark for Patient Safety Restraints and is owned by Robert Dan Allen.

TS7000dV naming corresponds to the following CE registered products: TruSystem 7000 (dV) - 1841049; TruSystem 7000 (dV) V - 1841083; Cable connection TS7000 (dV) - 1816914

The TS7000dV Surgical Table is intended to be used by clinicians and medically qualified personnel. This medical device is a regulated healthcare product which, pursuant to such regulation bears a CE mark. Hill-Rom recommends that you carefully read the detailed instructions for safe and proper use included in the documents accompanying the medical devices. The personnel of healthcare establishments are responsible for the proper use and maintenance of these medical devices. Hill-Rom reserves the right to make changes without notice in design, specifications and models. The only warranty Hill-Rom makes is the express written warranty extended on the sale or rental of its products. TRUMPF Medizin Systeme GmbH + Co. KG is a subsidiary of Hill-Rom Holdings, Inc.