

# Welch Allyn Insight™ Binocular Indirect Ophthalmoscope

Enhance fundoscopic exams at an affordable price.



- Insight is the lightest, wireless BIO on the market, weighing in at just over one pound
- Ergonomic Flip-Up™ optical system mount and on/off design eliminates the need for an overband and helps ease adjustment during an exam
- ColorSoft LED™ bright, white LED technology is optimized within the color spectrum and softer for enhanced patient comfort
- Rechargeable lithium battery pack power supply means no cords to impede your movement during an exam
- Unique light output control lever adjusts light intensity; an efficient alternative to the traditional solid state rheostat control
- Carries an exclusive 10-year industry-leading instrument warranty

# A full-featured device for enhanced workflow and efficiency, with the comfort and control you need.


Lightweight, rechargeable lithium battery pack provide up to 3.5 hours of on-time when fully charged

Headbands are ergonomically shaped, designed to match the contours of the head for maximized comfort

Padded, leather headbands are easy to adjust for improved stability

Power saving, Flip-Up optical system is on when pulled down/into operating position and powers off when flipped up

Light output control lever is integrated into optical system for convenient, hassle-free light adjustment

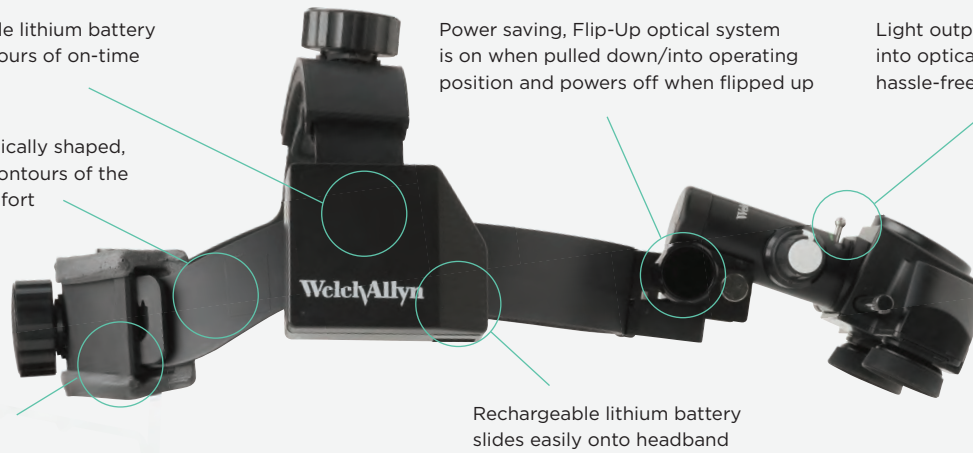

## Product Specifications:

Weight (device with rechargeable battery pack)	499 grams (1.1 lbs)
Rechargeable Battery Pack Life	3.5 hours when fully charged
LED Lamp Life	20,000 hours
Color Temperature	3100 °K
Filters	Neutral, Diffuser, Red-free, Cobalt Blue, Amber
Spot Sizes	Small, Intermediate, Large
Small Pupil Capability	Yes
IPD Range	48 to 74 m
Warranty	10-year instrument

## Ordering Information:

13500 Insight Binocular Indirect Ophthalmoscope with two Lithium Rechargeable Battery Packs, Charging Transformer and Wall Bracket

## Accessories:

- 75360 Lithium Ion Rechargeable Battery Pack and Charging Transformer
- 13550 Teaching Mirror with Case
- 13560 Wall Mount
- 13510 Charging Transformer
- 13570 Soft Carrying Case



# WelchAllyn®

Explore the entire Welch Allyn physical exam line at [welchallyn.com/exam](http://welchallyn.com/exam) or call us at **800.535.6663**.

Welch Allyn, Inc.  
4341 State Street Road, Skaneateles Falls, NY 13153 USA  
(p) 800.535.6663 (f) 315.685.3361

[WWW.WELCHALLYN.COM](http://WWW.WELCHALLYN.COM)

© 2015 WELCH ALLYN MCI2508 SM4150 REV A