



## Welch Allyn Stethoscope Specifications

### Welch Allyn Lightweight Double-Head Stethoscopes

When performance weight and economy are of primary concern, the Lightweight is the stethoscope to choose. At 4 oz, the Lightweight is one of the lightest stethoscopes available. It rests comfortably in your lab coat pocket or draped around your neck, and offers the superior acoustics that are the hallmark of Welch Allyn stethoscopes.

- General-use stethoscope
- Lightweight chestpieces with nonchill rims
- Adjustable chrome-plated binaurals
- Latex-free, flexible, single-channel tubing
- 2-year warranty

Lightweight Stethoscope  
Double-Head  
5079-73 Velvet Black  
5079-74 Poppy Red  
5079-75 Dawn Gray  
5079-76 Misty Blue



### Welch Allyn Uniscopes

Designed as a single-patient stethoscope for situations where cross-contamination issues are of concern, the Welch Allyn Uniscope disposable stethoscope enables you to make accurate diagnoses without compromising your primary stethoscope.

- General-use stethoscope
- For use in situations where universal precautions are mandated
- Single-head chestpiece
- Pediatric and anesthesia models available
- Latex-free tubing

Disposable Stethoscope  
17461 Adult, Yellow  
17462P Pediatric, Light Blue



PRODUCT DESCRIPTION	CARDIOLOGY GRADE		PROFESSIONAL GRADE		GENERAL PURPOSE	DISPOSABLE
	HARVEY DLX (28")	HARVEY ELITE DOUBLE HEAD (28")	PROFESSIONAL ADULT	PROFESSIONAL PEDIATRIC	LIGHTWEIGHT	UNISCOPE
PRODUCT NUMBERS	5079-321 Triple, Black 5079-322 Triple, Burgundy 5079-323 Triple, Navy 5079-324 Triple, Green 5079-325 Double, Black 5079-326 Double, Burgundy 5079-327 Double, Navy 5079-328 Double, Green	5079-125 Elite, Black 5079-270 Elite, Burgundy 5079-271 Elite, Navy 5079-284 Elite, Green 5079-125P Elite, Ped Blk	5079-135 Black 5079-137 Navy 5079-289 Blue 5079-139 Burgundy 5079-285 Forest	5079-145 Black 5079-147 Navy 5079-149 Burgundy 5079-287 Forest	5079-73 Velvet Black 5079-74 Poppy Red 5079-75 Dawn Gray 5079-76 Misty Blue	17461 Adult, Yellow 17462P Pediatric, Light Blue
PRODUCT TYPE	Double/Triple Head	Double Head	Double Head	Double Head	Double Head	Single Head
HEAD DIAMETER INCHES (mm)	1.9 (48.3)	1.9 (48.3)	1.75 (45)	1.2 (30)	1.9 (48.3)	1.9 (49)
OVERALL LENGTH INCHES (mm)	28 (711)	28 (711)	28 (710)	28 (710)	32 (812.5)	29.5 (750) 27 (686)
TUBING LENGTH INCHES (mm)	19.5 (495)	19.5 (495)	22.5 (571.5)	22.5 (571.5)	22 (558.8)	19 (483)
LATEX FREE (Y/N)	Yes	Yes	Yes	Yes	Yes	Yes
CHESTPIECE MATERIAL	Brass	Stainless Steel	Stainless Steel	Stainless Steel	Chrome Plated Zinc	ABS
DIAPHRAGM MATERIAL	Polycarbonate	Polycarbonate	Clear Polycarbonate Glass Epoxy Sheet	Iten Resiten G-10/FR4T	PVC	Polypropylene
DIAPHRAGM STYLE	Flat or Corrugated	Flat or Corrugated	Flat	Flat	Flat	Flat
ADJUSTABLE BINAURAL (Y/N)	Yes	Yes	Yes	Yes	Yes	No
EARTIPS (SOFT/HARD)	Soft/Hard	Soft/Hard	Soft/Hard	Soft/Hard	Hard	Hard
TYPE OF TUBING (SINGLE/ DUAL LUMEN)	Dual Lumen	Dual Lumen	Single Lumen	Single Lumen	Single Lumen	Single Lumen
ANTI-CHILL RING (Y/N)	Yes (Brass)	Yes	Yes	Yes	Yes	No
SHIPS ASSEMBLED (Y/N)	Yes	Yes	Yes	Yes	No	Yes
PACKAGING	Printed Box	Printed Box	Printed Box	Printed Box	Plain White Box 1-piece	Plastic Bag

**WelchAllyn**  
Advancing Frontline Care™

4341 State Street Road, PO Box 220, Skaneateles Falls, NY 13153-0220 USA  
800.535.6663 www.welchallyn.com

© 2012 Welch Allyn SM2372 Rev E

LIGHTWEIGHT

DISPOSABLE


# Welch Allyn Stethoscopes



**WelchAllyn**  
Advancing Frontline Care™

## Welch Allyn Stethoscopes *Hear the Difference*

Developed during a 30-year association with W. Proctor Harvey, MD, recognized “Master” by the American College of Cardiology and one of the world’s foremost experts in auscultation, Welch Allyn stethoscopes bear the same reputation as all of our diagnostic instruments—rich in innovative features and design, unsurpassed in exceptional quality and performance, and representative of a long and collaborative association with the medical profession.

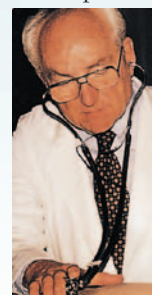

Welch Allyn stethoscopes embody a proven history of success. And continual research and refinements will ensure they remain the choice of healthcare professionals in the future.

Welch Allyn products are available exclusively through authorized dealers. For help locating a dealer or product, contact a Welch Allyn representative at 800.535.6663, or visit [www.welchallyn.com](http://www.welchallyn.com).



### Welch Allyn Harvey DLX™ Stethoscope

Both the double- and triple-head Harvey DLX stethoscopes are specifically designed to accentuate important patterns and to enable you to hear cardiac, pulmonary, and vascular nuances without relying on difficult “variable pressure” techniques. The Harvey DLX is equipped with dual-lumen tubing that creates two uninterrupted sound channels and works with the optimally weighted brass chestpieces to transmit a full range of physical sounds.

- Full-range, cardiology/pulmonology stethoscope
- Available with double-head (brass bell and flat diaphragm) or rotatable triple-head (brass bell and flat or corrugated diaphragm) chestpieces
- Pediatric chestpieces available
- Chrome-plated brass chestpieces and binaurals
- Rotatable binaurals and interchangeable comfort sealing eartips
- Triple-leaf binaural spring encased in polyurethane allows size adjustment while preventing breakage
- Latex-free
- 10-year warranty

Harvey DLX Double-Head  
5079-325 Black  
5079-326 Burgundy  
5079-327 Navy  
5079-328 Forest Green  
5079-325P Black, Pediatric

Harvey DLX Triple-Head  
5079-321 Black  
5079-322 Burgundy  
5079-323 Navy  
5079-324 Forest Green



### Welch Allyn Harvey™ Elite® Stethoscope

The Harvey Elite offers updated styling along with the superior comfort and auscultation capabilities that cardiologists have come to rely on. It is equipped with dual-lumen, latex-free tubing that creates two uninterrupted sound channels and works with the optimally weighted, rugged stainless steel chestpiece. The Elite’s diaphragm is designed to accentuate important rhythms, enabling you to better hear cardiac, pulmonary, and vascular nuances.

- Full-range cardiology/pulmonology stethoscope
- Double-head (bell and flat diaphragm) stainless steel chestpiece with nonchill rim
- Pediatric chestpieces available
- Rotatable stainless steel binaurals and interchangeable comfort sealing eartips
- Triple-leaf binaural spring encased in polyurethane allows size adjustment while preventing breakage
- Choice of tubing length
- Latex-free
- 10-year warranty

The Harvey Elite (22")  
5079-122 Black

The Harvey Elite (28")  
5079-125 Black  
5079-270 Burgundy  
5079-271 Navy  
5079-284 Forest Green  
5079-125P Black, Pediatric



### Welch Allyn Professional Stethoscope

With minimal weight, custom-fitting binaurals, and eartips, the Professional is one of the most comfortable and versatile general-use stethoscopes manufactured today. Its acoustics allow for accurate auscultation of cardiac, pulmonary, and vascular sounds. And the chestpieces in the adult, pediatric and infant versions of the Professional enable you to select a stethoscope that best meets your specific needs.

- General-use stethoscope
- Stainless steel double-head (bell and flat diaphragm) chestpiece with nonchill rim
- Available in pediatric and adult models
- Rotatable binaurals and interchangeable comfort sealing eartips
- Latex-free
- 5-year warranty

The Professional  
5079-135 Black  
5079-137 Navy  
5079-139 Burgundy  
5079-285 Forest Green



# HARVEY™ DLX

# HARVEY™ ELITE®

# PROFESSIONAL